APP Driven Technology powered by

Prof. Dr. Kafka

USER MANUAL

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FEEL THE CHANGE

1.	Warning & Safety Information	01
2.	Disclaimer & Terms of Use	03
3.	Intended Use	04
4.	Congratulations ...	05
5.	Prof. Dr. Wolf A. Kafka	06
6.	Scope of Delivery	07
7.	Usage of the KLOUD+ Application	08
8.	Initial Commissioning	09
9.	Control Unit	10
10.	CENTROPIX App	15
11.	Application	17
11.1.	Full-body Application	18
11.2.	Local Application	19
11.3.	Examples of Use	21
11.4.	Programs for Application	22
12.	SCAN-X	23
13.	Battery	23
14.	Repair and Cleaning	25
15.	Disposal	27
16.	Troubleshooting	28
17.	Specifications	29
18.	Warranty	31

1. Warning & Safety Information

Be sure to observe every detail of this User Manual and always keep it close to your **KLOUD+** applicators. Before use, this User Manual must be read and understood. Of course, this assumes that you are able to read. Laws require a reading ability of at least 5 years. Otherwise, the operation must be carried out by a suitable person.

These operating instructions apply exclusively to **KLOUD+** applicators, and contain important instructions for commissioning and handling. Please also contact your local CIA (CENTROPIX Independent Affiliate) for a brief introduction to the **KLOUD+**.

For translations of the content into a foreign language, which is not available in this user manual, the German language version is the binding original.

-  Children shall not play with the appliance.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- KLOUD+** applicators do not need to be maintained or calibrated, however, it is always advisable to keep the **KLOUD+** in order at all times to avoid possible accident hazards! In case of damage, these parts are to be replaced only by the manufacturer (see warranty service).
- The warranty is only guaranteed when using the **KLOUD+** up to a maximum body weight of 120 kg (265 pounds) and when used on a flat, hard surface. Please place the **KLOUD+** on your body if these conditions cannot be met.

- Protect yourself and your **KLOUD+** applicators from damage by using them only on stable, non-slip and smooth surfaces. The **KLOUD+** applicators are made of very high-quality materials, but sharp objects can damage the material.
- When unfolding and folding the **KLOUD+**, make sure that you do not damage the upper part of the **KLOUD+**.
- Make sure that the head part of the applicators (with the control unit) does not bounce against hard objects and is not subjected to heavy loads (e.g. do not step on the head part or place heavy objects on it).
- Always place the applicators in the appropriate pocket with the head section facing upwards. This prevents the control unit from being damaged when the bag is put down.
- Make sure that the mains voltage is identical to the specification on the power supply unit before connecting the **KLOUD+** to the mains for charging. Only use the enclosed power supply unit, adapter and charging cable.
- Since the **KLOUD+** is battery powered, do not dispose of it in household waste. Please observe the RoHS standards listed in this User Manual.
- The **KLOUD+** must not be exposed to any heat source and must not be thrown into a fire.
- The **KLOUD+** is not waterproof.
- Applicable regulations and laws must be observed by professional users (e.g. doctors, physiotherapists, etc.).
- The **KLOUD+** has undergone and passed rigorous safety testing. CE and FCC standards ensure safe use.

2. Disclaimer & Terms of Use

⚠ As a manufacturer, we follow the obligation not to make any healing statements about the **KLOUD+**. The statements and information made about the effects of the device are source-based, scientific, reputable and trustworthy studies, as well as experience reports from users. We dissociate ourselves from any healing statements and promises of healing, as well as from statements of medical efficacy within the meaning of the Therapeutic Products Act (TDA) and the Act against Unfair Competition (UWG).

Thousands of users have confirmed the positive properties of **KLOUD+** and no longer want to do without their **KLOUD+**. Negative effects are not yet known. Despite all the positive empirical research and results, as a wellness and fitness company, we advise you to consult a licensed physician for health issues or concerns.

The use of the device has not been evaluated by government regulatory authorities (e.g. FDA, etc.). The descriptions are not intended to diagnose, treat, cure or prevent any disease. No claim or opinion about general health is intended as medical advice and should not be construed as such.

Please consult a licensed physician before starting any exercise, nutrition or other wellness and fitness program.

CENTROPIX Global AG and its subsidiaries, independent partners and sales partners are not liable for the use of **KLOUD+**. USERS take full responsibility for the health, wellness and fitness decisions they make and pursue.

⚠ The following users may only use the **KLOUD+** applicators after consultation with their doctor:

- If you are suffering from a chronic degenerative diseases or have recently undergone surgery
- People with a medication regime and constant medication intake.
- People who wear a pacemaker, electronic implants and subcutaneous medical devices.
- If you are a pregnant women.
- People who have undergone foreign organ or cell transplants and have therefore been prescribed immunosuppressants by their doctor.

3. Intended Use

⚠ The **KLOUD+** are not medical products; they are wellness & fitness products that may support your general performance, recovery, wellbeing and quality of life by means of clearly defined, low-frequency pulsating electromagnetic field conditions.

- The operating instructions are for information and handling purposes only and are not to be considered as medical advice.
- If you have a medical problem, we recommend that you consult a doctor and follow the previous warnings and disclaimers.

CENTROPIX **KLOUD+** products may not be used for any purpose other than those specified in these operating instructions and their intended use.

Any other use or application of the **KLOUD+** devices voids the guarantee. Furthermore, the manufacturer assumes no liability for any resulting damage.

4. Congratulations ...

... to your pulsating electro-magnetic field system **KLOUD+**!

By purchasing this product, you are doing a lot of good for your health. This is because the **KLOUD+** may support your general performance, recovery, wellbeing and enhance your general quality of life by means of electromagnetic impulses and can thus support your self-healing powers.

What motivated us to develop the **KLOUD+** is today's difficult time, which is very difficult for us all. The general performance decreases, the joy of life decreases and the quality of life is somehow lost. Most people are not ill, but in many cases overworked and burnt out. The body is crying out for help and just withdrawing does not improve the situation as we all know. The solution can be to help the body to help itself. To bring life back into the balancing center - completely under the motto: Be the change.

In collaboration with the physicist and neurophysiologist Prof. Dr. W. Kafka, we have created this unique product. The **KLOUD+** was developed as a **mini** and **maxi** version and these two applicators differ only in size. Both variants therefore have advantages in terms of mobility and its areas of application. By the way: This user manual is therefore equally valid for the **KLOUD+ mini** as well as for the **KLOUD+ maxi**. These two applicators support your health with a very wide range of applications, which we hope you enjoy.

Your CENTROPIX team

5. Prof. Dr. Wolf A. Kafka - the mastermind of PEMF technology

Prof. Dr. Wolf A. Kafka, the inventor of the unique, patented signal configurations of the **KLOUD+** Technology is a theoretical physicist and neurophysiologist who has been dealing with the subject of electromagnetic field therapy for more than 30 years. And, by many recognized as the Father of PEMF technology.

He has been accredited with numerous technical developments, patents and award-winning documentations that garner his claimed achievements and accomplishments.

In 1998, Prof. Dr. Kafka patented a physical method which application is aimed at supporting natural self-healing powers, i.e. the human body «inner doctor». This was achieved by constructing special electromagnetic pulses.

Thousands of users reported marked improvements in their general quality of life, better sleep, and much more with these previous devices. And it is no surprise that it became a new standard in wellness and fitness applications worldwide.

In 2019, Prof. Dr. Wolf A. Kafka took into consideration all these experiential users' feedback to build brand-new, improve and broader ground signal including other technical advantages which led to a brand-new 2021 international patent.

With all these years of research, development and successes Prof. Dr. Wolf A. Kafka once again surprises the world by introducing a new improved form of P.E.M.F. called P.E.M.A (= Pulse Electromagnetic Amplitude). This latest invention is what you have in your hands: the CENTROPIX **KLOUD+** system. It does not only set a new standard in customer benefits and improvement, also sets new trends and standards in user friendliness, user experience, design and quality.

6. Package Contents

		KLOUD+ complete	KLOUD+ maxi	KLOUD+ mini
KLOUD+ Suitcase		✓	✗	✗
KLOUD+ Designbox In 2 sizes depending on the product		✗	LARGE	SMALL
KLOUD+ maxi Applicator for whole-body application		✓	✓	✗
KLOUD+ mini Applicator for partial body application		✓	✗	✓
SCAN-X To make the signal audible		✓	✓	✓
Power supply & adapter + charging cable		2 SETS	1 SET	1 SET
Elastic band To fix the applicators		2 PIECES	1 PIECE	1 PIECE
CENTROPIX App For smartphone or tablet		DOWNLOAD	DOWNLOAD	DOWNLOAD
Bag		LARGE & SMALL	LARGE	SMALL
Lint Glove		✓	✗	✗
User Manual		✓	✓	✓
Quick Start Guide		✓	✓	✓

7. Usage of the KLOUD+ Application

The definition of physiology is as follows: A branch of biology concerned with the normal functions of living organisms and their parts. Physiological effects of electromagnetic fields depend on various factors:

- The variation in intensity over time
- From the spatial distribution
- Frequency bandwidths
- Amplitude curves
- Clarity of the magnetic fields

The 8-coil applicators from KLOUD+ take this into account. The variation offers flexible application options adapted to individual body applications with different depth effects, with the simplest handling.

8. Initial Set-Up

1

Remove the **KLOUD+** devices from the design box and place them on a flat surface.

2

Before using the applicator for the first time, charge the batteries. To do this, connect the power supply unit to the respective applicator with the USB cable (see chapter Battery). Make sure that you first click the power plug that corresponds to the power supply in your country into the adapter. You can use the other adapter in other countries, e.g. when travelling.

The power supply unit and plug must be in undamaged original condition. Contact CENTROPIX for replacement of damaged parts immediately!

Info: When removing the magnetic charging cable, the dongle can come off. Disconnect the dongle from the charging cable and plug it back into the socket.


3

The charging time is approx. 4 hours. The battery is full when the blue L.E.D. the charging symbol goes out at the battery, and the green L.E.D. on the charging cable lights up.

Press the «Power» button above the keyboard. The «Power» L.E.D. Lights up immediately.

After approx. 15 seconds the device is ready for operation and the «Ready» LED lights up. The blue LED and the green P1 LED light up.

9. Control Unit

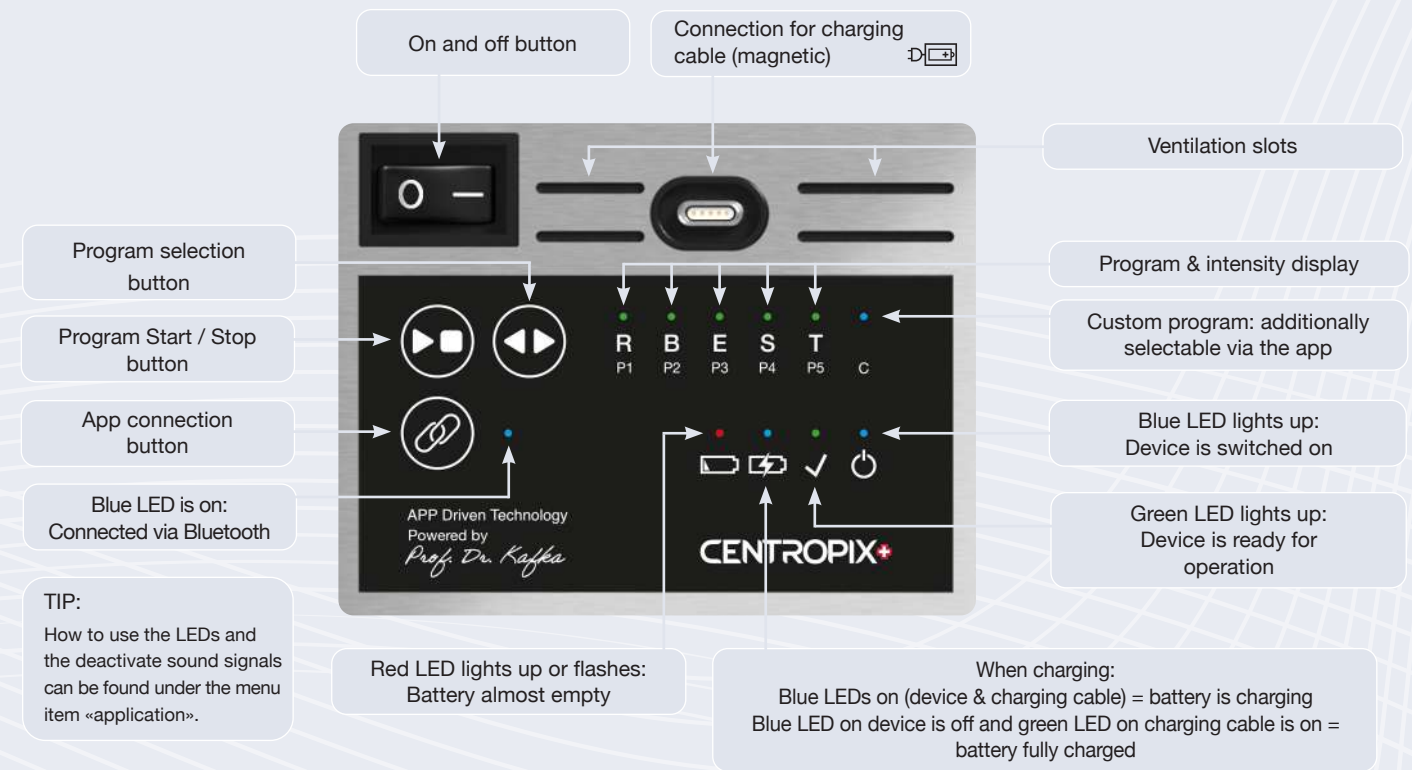The control unit is the central power station for your **KLOUD+**. You can conveniently transport the complete applicator including the control unit and use it at different locations.

Make sure that the control unit (headboard) is not exposed to strong vibrations and stress. Do not stand on the head unit!

Next to the keyboard, there are ventilation slots that serve to cool the computer unit. Never cover these slots!



The applicators are switched on and off by a mechanical push button. 


We recommend using the CENTROPIX App to control the applicators to maintain all functions and the full comfort of the **KLOUD+** applicators.

Privacy

The use of the devices and the app are subject to the strict data protection laws of the EU/USA. No personal data will be collected, recorded, or passed on.

Regulation of the magnetic field intensity

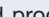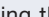Start the desired program. Check the set intensity by pressing and holding the start-stop button . After approx. 2 seconds the LEDs from P1 to P5 indicate the intensity. Continue to hold the start-stop key and use the Program Select key  to gradually select the desired intensity:

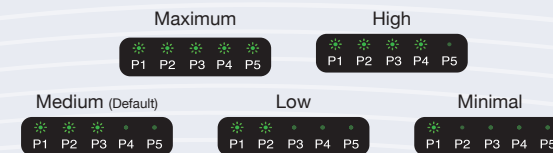Thus, the intensity applies to all selected programs until it is changed again (remains even after the device is switched off).

After releasing the start-stop button , the display is back in standard mode and the selected program continues.

This control function is, of course, also possible via the app.

Character Explanations

-  FCC-compliant USA
-  CE-compliant EU
-  Note
-  User manual
-  Manufacturer
-  Do not iron
-  Do not wash
-  Do not tumble dry
-  Do not dry clean
-  Disposal not in household waste
-  Recyclable
-  RoHS compliant
-  Battery charging
-  Battery empty
-  Caution / Attention
-  Correct
-  Incorrect
-  Direct current

Standby mode

- When the **KLOUD+** is started, the blue LED flashes  for approx. 15 seconds and the **KLOUD+** is ready to be used. If no further button is pressed, the »standby mode« automatically starts after 5 minutes.
- The »standby mode« requires very little battery. Therefore, with a fully charged **KLOUD+**, the battery will last up to 3 weeks in »standby mode«.
- If the »on-off button«  is turned off, the device is NOT in »standby mode«.
- If the »start-stop button«  is pressed for 3 seconds, the device switches to »standby mode«. The blue LED flashes  during shutdown and then switches off.
- To activate the device from »standby mode«, press the »start-stop key«  for 1 second and the device is immediately ready for use.
- Since all LEDs are switched off during »standby mode«, it is not possible to see whether there is still enough battery for an application. To check the battery status, switch the »on-off switch«  off and then on again. If the blue LED  flashes, the device boots up and is ready for use again after 15 seconds. Otherwise please charge the battery.
- During »standby mode«, you cannot connect the app to the applicator.

Key Combinations at a Glance

Start program

-  Set the on-off switch to „I“
-  Select the desired program with the program selection button
-  Press the start-stop button once

Magnetic field intensity

-  Start program
-  Hold the start-stop button for 2 seconds to see the intensity level
-  While holding the start-stop button, select the desired intensity with the program selection button.

Standby mode

-  Hold the start-stop button for three seconds
- Or: Do not press any button for 5 minutes (Exception: program is running)

Do-not-disturb Mode

-  Hold the program selection button
-  and the start-stop button simultaneously.
- Now, only the “ON” LED is on.

Check battery level

-  Hold the program selection button for 3 seconds

Activate Bluetooth

-  Press the Bluetooth button once

Deactivate WiFi in the app

Open the App on your smartphone > Menu > Settings > WiFi configurations > Turn on the WiFi toggle key > Enter WiFi name and password > Connect > Pair.

10. CENTROPIX App

CENTROPIX App Download

You will find the CENTROPIX app for your mobile phone as a free download at www.centropix.com/app as well as in the App Store or the Google Play Store.

You are therefore comfortable and numerous operating options are available.

Updates

KLOUD+ applicators are sustainable and therefore designed so that product developments can be transferred via the app with an appropriate software update. Updates can take between one minute, and in some cases up to 4 hours.

To perform an update, you must first connect your applicator and set up the WiFi function in your app: Menu > Settings > WLAN Configuration > Slide WiFi on > Enter your WiFi name and your password > Connect > Pair. Important: WiFi will only be utilized by the **KLOUD+** to perform updates.

You are informed whether an update is available in Settings > Updates or through a button highlighted in blue on the dashboard. This button only appears if the WiFi function is set up.

To prevent an update from being aborted due to the battery status, we recommend that you first fully charge the Applicator. In order to successfully install the update, please ensure that your **KLOUD+** has been charged at least 40%.

The blue «Connect LED»  flashes during the update. When it lights up continuously, the software is installed on your Applicator. The update process must not be interrupted! To complete the update, you must restart your **KLOUD+**. Your Applicator is now up to date.

You are also welcome to run your application during an update.

Control of multiple devices

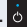
Multiple applicators can be operated via a smartphone, tablet, etc. by using the CENTROPIX App.

11. Application

The **KLOUD+** **maxi** is suitable for whole-body applications, whereby, like the **KLOUD+** **mini**, it can also be ideally used for local applications when folded.

Tips for your optimized **KLOUD+** application:

- Be sure to consume a glass of natural water before and after use, but do not consume coffee or tea before use.
- Do not consume tobacco after using **KLOUD+**.
- Make yourself comfortable: Take a comfortable position in comfortable clothes.
- Ensure the right environment: choose a place with a comfortable temperature without any stress or distraction, and use pleasant relaxation music if possible.

- 1.) Switch on the desired applicator with the power button.  The blue power LED lights flashes. 
- 2.) Wait until the OK LED lights up. 
- 3.) Select the desired program using the program selection button  (P1-P5 by using the keyboard or others by using the Centropix operating app).
- 4.) Press the start button  to start the selected program. The programs duration is normally 15 minutes. While a program is running, the corresponding program number LED flashes.
- 5.) After the program has finished, the program LED stops flashing and lights up permanently again when another program is selected. An acoustic signal sounds at the end of the program.

«DO NOT DISTURB MODE»

To switch off the LEDs and the acoustic signal during use, i.e. to activate the «do not disturb mode», please press the two buttons «Select»  and «Play» . In a dark environment the lighting can be reduced to one LED («ON») (this cannot be switched off). By pressing any key you exit the «do not disturb mode» and the LEDs will be displayed again.

11.1. Full-body Application

The general full-body application is carried out by default with the **KLOUD+** **maxi** application module with an average intensity of 100 μT . When the applicator is folded three times, the maximum average intensity is 180 - 200 μT . For comparison: The natural earth magnetic field strength is approx. 35 - 70 μT depending on the latitude.

The field distribution extends slightly beyond the applicator surface. It is therefore of little importance whether head or feet extend beyond the ends of the applicator during usage.

Also, it does not matter whether the top or bottom of the applicator is used since the signal is generated identically on both sides of the coil.

The arrangement of the coil elements is symmetrical. Therefore, the applicator can be placed, e.g. if you want to place it in a bed, with the control element (foam element) at the head end as well as at the foot end.

People with limited physical movement can be «covered» with the **KLOUD+** **maxi** both in a sitting and lying position. This way you can enjoy the application on the sofa, bed, etc. comfortably and conveniently.

We recommend two Full body applications daily (one in the morning one in the evening).

However, do not use more than one **KLOUD+** applicator at the same time.

11.2. Local Application

When unfolded, the **KLOUD+** mini has been designed to fit the standard size of a normal office chair. Thus, the **KLOUD+** mini finds its place in tangible life without major disruptive factors, because it is extremely handy and mobile thanks to its foldable design. It can be integrated into everyday life at the workplace, in the car, on a plane, and much more, without taking up additional time. For example, there is the option to use the **KLOUD+** mini applicator while sitting on a chair or to cover oneself with this smaller applicator.

When folded, the **KLOUD+** mini allows more intensive applications.

The average intensity (microtesla) of both applicators is about 100 μ T. This varies within the programs. Depending on the folding of the applicators, the average intensity increases. Thus, with 3-fold folding, between 180 - 200 microtesla are achieved.

When folding, make sure that the squares lie directly on top of each other.

The **KLOUD+** mini is simply held at/on the body part to be treated. It can also be easily fixed to the back of a chair with the included rubber band.


The large applicator, the **KLOUD+** maxi, can also be folded like the **KLOUD+** mini with the same effects.

Choose the appropriate program (description P1-P5, see «Programs for application»; P6 and others Programs: Operation and information via the app) for your current situation and the desired effect.

Application of the **KLOUD+** is possible several times a day.


IMPORTANT! If the folding is not done correctly, it is possible that the intensity of the magnetic field will be reduced or neutralized. The intensity of the magnetic field can only increase if the KLOUD is folded correctly. The control unit is always positioned at the top left when folded correctly.

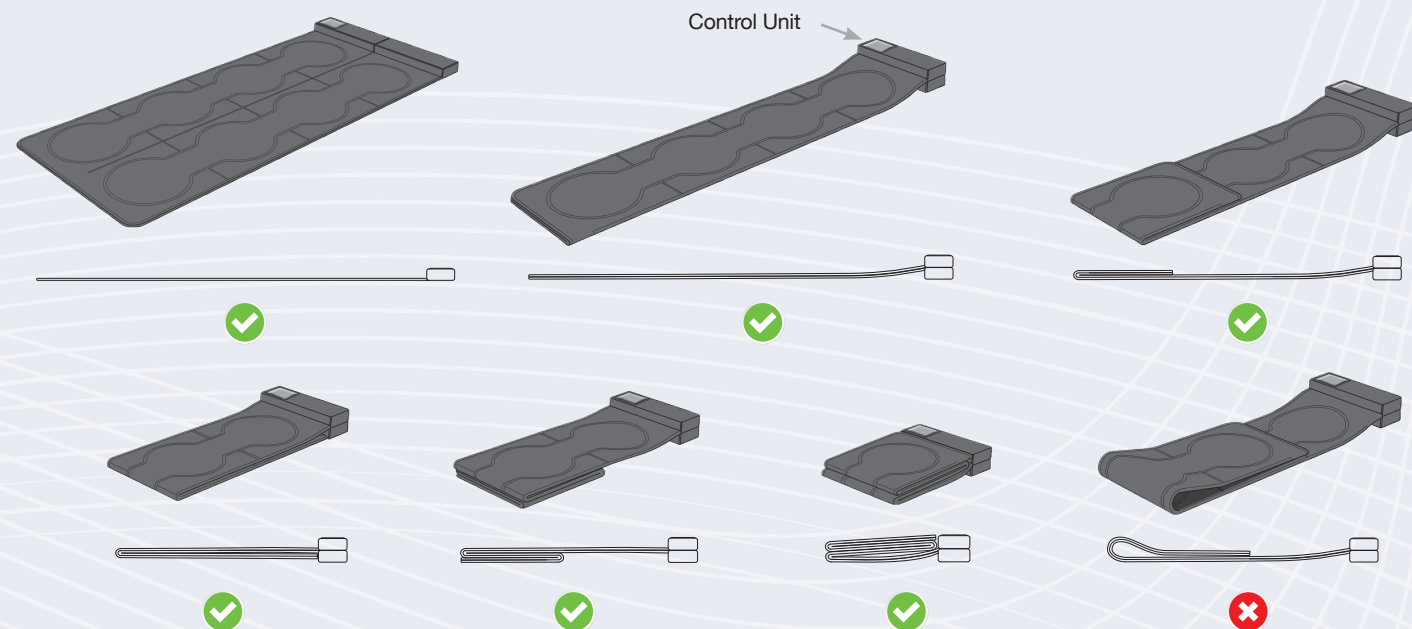
11.3. Examples of Use

KLOUD+ maxi

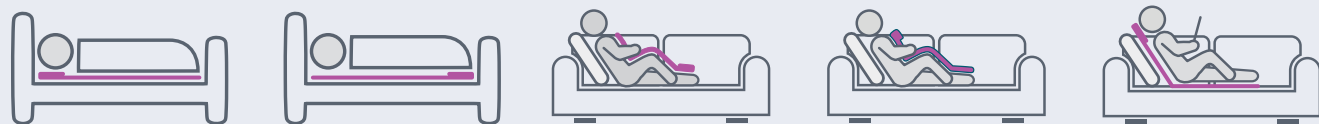KLOUD+ mini

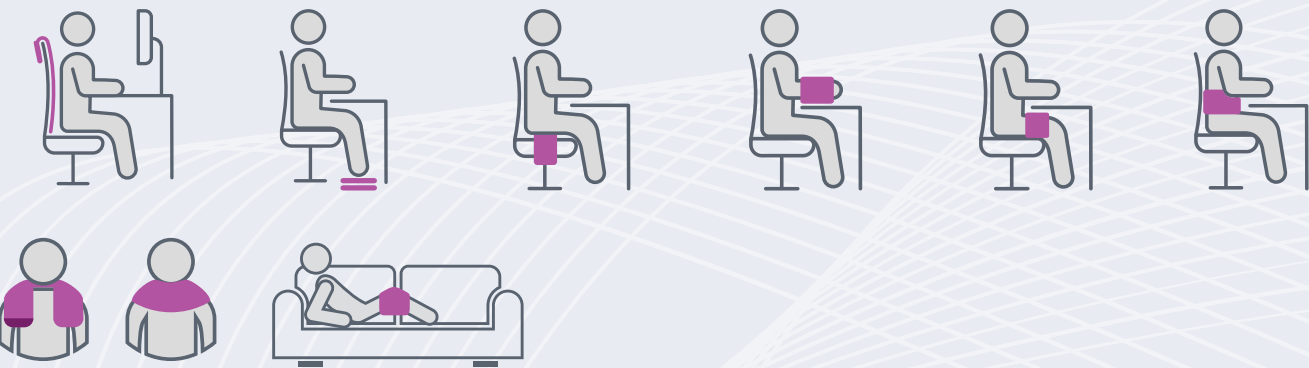11.4. Programs for Application

The KLOUD+ devices contains five fixed programs with different signal configurations:

P1

RELAX

Relaxation and entrainment of body, mind and soul. Put yourself in a relaxed state and unwind. Promote your well-being and find inner peace, and optimize a restful sleep.

P2

BALANCE

Restore your natural balance and get into a good flow. Balance means feeling good and thinking positively with confidence, peace and serenity. Gain peace and quality of life.

P3

ENERGIZE

A boost when you want more life energy! Experience pure "KLOUD+ power", feel alive, awake, motivated, and strengthened. Entrain and activate vitality from within. Take advantage of your bodies natural power! Rhythm like the heartbeat and breath rate are kicking your vitality into a higher level.

P4

SYNCHRONIZE

The ideal introduction to the use of the KLOUD+. Synchronize your body functions. Use this program in combination with relaxation habits and rituals to improve mental clarity, emotional harmony and synchronize your heartbeat.

P5

TRANSFORM

Help the body help itself. Keep your vital forces turned up, and experience a positive transformation - push your boundaries!

* The KLOUD is not intended for the detection, prevention, monitoring, treatment or alleviation of disease. In particular, the KLOUD is not a medical device, nor are any other medical services provided through the KLOUD. The applications provided through the KLOUD should not be relied upon to make medical decisions or to diagnose or treat any disease. If you have health problems or questions about your health, always consult your physician. The programs are general recommendations based on experience. Individual experiences and results may vary. All KLOUD programs described, and their properties are based on empirical evidence that has not been evaluated by a competent authority - such as the U.S. FDA (Food & Drug Administration).

More programs are available through CENTROPIX App.

12. SCAN-X

With the SCAN-X module, you can display the electromagnetic field by means of an acoustic signal and is therefore suitable as an effective demonstration for those interested.

To do this, switch on the Scan Module by pressing the button, keep it pressed and hold it in the area to be tested (L.E.D. lights up).

The more intensive the signal to be measured, the louder the sound of the SCAN-X.

It is best to use the program P3: «Energization».

If the battery is low, the red L.E.D. is only displayed during the first 2 seconds.

Suppose the battery capacity is less than 10%, the red L.E.D. lights up.

Suppose the battery capacity is less than 5%, the red L.E.D. flashes.

Changing the battery of the SCAN-X

Open up the SCAN-X Case by pushing apart the upper and lower case (i.e. by using a screwdriver). Replace the CR-2032 battery with a new model.

13. Battery

Do not use the KLOUD+ while it is charging. Both KLOUD+ applicators are equipped with a high-quality lithium-ion battery. A new, fully charged battery lasts for over 4 hours of use or approx. 17 applications of 15 minutes each. Depending on the battery care or age of the battery, the number of applications per battery charge may be lower.

Check battery status

Press the program selection button for 3 seconds. Depending on the battery status, the LEDs P1 to P5 light up (example: All LEDs from P1 to P5 light up means that the battery is full).

The charging time is about 4 hours. During the charging process, the charging L.E.D. «Charging battery» is lit.

Exception: If the battery is very low, it can take one to two minutes to start the device, then the blue power L.E.D. Lights up, but the system cannot start. In this case, please connect the charging cable.

As soon as the battery is charged, the charging process is automatically terminated and the blue charging L.E.D. goes out.

When the battery status is approx. 20%, the L.E.D. «Low battery» lights up. Then approx. 2 more therapies can be performed. At a battery level of approx. 10%, this L.E.D. starts flashing. You can then perform approx. 1 therapy before the battery needs to be recharged.

The 1.8 m long charging cable and the magnetic connections of the KLOUD+ maxi and the KLOUD+ mini ensure convenient charging. Connect the charging cable to the applicator and the green L.E.D. on the charging cable changes to blue. If the L.E.D. does not change to blue, the battery is already fully charged. In addition, the blue charging L.E.D. on the applicator keyboard lights up during the charging process.

Fold the KLOUD+ and open the leather flap on the side of the headboard. Fold the flap underneath and carefully pull out the white plug. Disconnect it by pressing the release clip. Pull out the battery which is the short part of the plug.

Only the original battery and the supplied detachable power supply unit with the original charging cable may be used. If another battery or charging cable is used, the device can be damaged, e.g., the warranty becomes void.

Replacement batteries or charging cables can be ordered through your CENTROPIX consultant or online on the CENTROPIX Partner Page.

14. Repair & Cleaning

Repair

Before you send your CENTROPIX product for repairs to CENTROPIX or to an authorized specialist company designated for this purpose, you must first carry out the general troubleshooting (see chapter „Troubleshooting“). For this purpose, we urge you to contact your supervising CENTROPIX partner or CENTROPIX by e-mail at rma@centropix.com.

It is important to us to be able to provide optimal service even in the event of a repair. Therefore, in the event of a repair, a so-called „First Troubleshooting“ must be carried out by e-mail under rma@centropix.com.

General rules:

- Replace damaged parts immediately, which you can get from CENTROPIX.
- Contact the appropriate service point to initiate the repair process.
- To ensure that you can safely ship your **KLOUD+** for repair, we recommend that you ship it in its original box.

Cleaning

Regular cleaning is recommended. Make sure that the device is unplugged during cleaning.

- Cleaning and user maintenance shall not be made by children without supervision.
- Clean and disinfect the **KLOUD+** thoroughly before it is used by another person.
- Remove lint, hair and small particles from the applicator with the included cleaning glove. Brush against the line to catch the lint, with the stroke remove the lint from the glove.
- Abrasives, acidic or alkaline cleaning agents can damage the **KLOUD+** and must therefore not be used for cleaning.
- Sharp-edged and pointed objects can damage the **KLOUD+** and should therefore not be used during cleaning.
- When cleaning and disinfecting the **KLOUD+**, make sure that no liquids enter the control unit. The **KLOUD+** applicators and associated components are neither water-proof nor splash-proof! The pads are not washable.

BE THE CHANGE

15. Disposal

This CENTROPIX product, including all parts contained therein (e.B cables, wires, etc.), comply with the requirements of Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment (RoHS).

When disposing of the appliances or batteries, observe local regulations. Batteries do not belong in the household waste and must not be thrown into the open fire, otherwise there is a risk of explosion. Be sure to ensure that the batteries are disposed of properly.

If you need further information, please contact your advisor.

16. Troubleshooting

In case of an error message or suspected malfunction, check the following points first. Errors are generally indicated by acoustic beeps and two flashing LEDs.

ERROR	POSSIBLE CAUSE	REMEDY
No L.E.D.s lighting up after 15 seconds boottime (after turn on power).	Battery empty	Charge battery with charging cable.
Power L.E.D. and indicator L.E.D. flash and a signal tone sounds.	<ul style="list-style-type: none"> • Computer error • The cable connection in the device interrupted 	Switching off and on again.
After connecting to the app (blue L.E.D. on) the L.E.D. flashes and the connection is interrupted. (Message in the CENTROPIX-App).	Bluetooth connection error between app and KLOUD+ .	Press the Connect button again, if necessary, confirm the new connection on the smartphone/tablet.
Requested program not available (custom programs).	Program is expired. Program not booked.	Book the desired program (again) (via App) or use one of the 5 pre-installed programs.
The device is switched on but does not start the program.		Establish a connection between applicator(s), update starts automatically.
The green «OK-LED» flashes.	Contact defect in the applicator.	Device must be sent to CENTROPIX.

If an error cannot be corrected or can be specifically located, please proceed as described in the chapter „Repairs“.

17. Specifications

Category	KLOUD+ maxi	KLOUD+ mini
Item No.	1128	1129
Measurement (LxW)	164 x 67 cm (64.57" x 26.38")	98 x 43 cm (38.58" x 16.93")
Weight	2,2 kg (4.85 lb)	0,99 kg (2.18 lb)
Serial No.	Fold up the KLOUD+. There is a leather flap on the side of the head section that leads to the battery compartment. Open this leather flap and fold down the lower flap underneath. There you will find the serial number of the KLOUD+. You can also find the serial number on the last page of this User Manual.	
Operating voltage	12-15, 1VDC	
Surface material *)	Front side 100% polyester, back side: 100% PVC PES, Flexathen FT/SPE 33 FR, polyethylene foam Flame retardant, stain resistant FMVSS 302 / ISO 75200 / BS 5852 Source 1, EN 1021-1&2	
Power consumption	250 - 650 mAh (when charging the battery at up to 1,7A)	
Coils	8	
Power consumption	Max. 8,5 W	
Power output (on the applicator)	Max. 6 W	
Maximum flux density of the application	Mat laid out: ca. 100 µT 3 x folded to: ca. 200 µT	
Connections	Magnetic charging socket	

*) The surface material can change depending on the batch, but this does not affect the high quality level. The colors of the fabrics can also vary slightly.

Category	CLEANING GLOVE
Measurement (LxW)	approx. 21 x 14 cm (8.3" x 5.5")
Weight	80 g (2.82 oz)
Material	60 % Polyamide 40 % Cotton

Category	CHARGER
Item number	1131
Measurement	75 x 40 x 30 mm (3.0" x 1.6" x 1.2")
Weight	140 g (incl. cabel) (4.94 oz)
Model	FOX12 USB
Adapter	EU (10 g) / US (0.35 oz)
Input	100-240 VAC ~ 50/60 Hz/300-150mA
Output	5 VAC - 2.2 A
Part-No.	5317-FW8000/USB
Spec.-No.	15.240VAC
Use	Indoor only
Material	PC / ABS V0 125° C / 257 F°

Category	AKKU
Capacity	4000 mAh
Minimum capacity	3900 mAh
Nominal voltage	3,6 V - 3,7 V
Final charge voltage	4,2V ± 0,05
Final discharge voltage	2,5 V
Max. discharge current (constant)	45A with switch-off at 80° C / 176 F°
Max. discharge current (pulse)	35A without shutdown
Max. discharge current (pulse)	approx. 70 A
Diameter	21,15 mm ± 0,15 mm (0.8327" ± 0.0059")
Height	70,15 mm ± 0,2 mm (2.7618" ± 0.0079")
Weight	70 g ± 1 g (2.47 oz ± 0.04 oz)
Charging procedure	CC-CV

Category	USB-KABEL
Item number	1116
USB type	USB Typ A plug USB 3.1 Typ C plug Magnetic connector & dongle
Jacket Material	Soft - PVC
Length	190 cm (74,80")

Category	SCAN-X
Item number	1130
Weight	60 g (2.12 oz)
Measurement	approx. 70 x 40 x 13 mm (2.76" x 1.58" x 0.5")
Operating voltage	3 Volt
Power supply	Button cell CR 2032
Measuring range	10 Hz - 1 KHz
Signal	Sound

18. Warranty

Warranty period

CENTROPIX grants users a warranty of 3 years (36 months) from the date of purchase for the **KLOUD+** applicators. LEDs and battery have a warranty of 6 months.

Warranty extension

At the time of purchase, the warranty can be extended to 60 months (5 years) at an extra charge. A subsequent warranty extension is not possible.

Prerequisite for warranty services

A warranty is only guaranteed to first-time users who have purchased the **KLOUD+** applicators via the official store of CENTROPIX. This warranty covers material or manufacturing defects only. In case of non-intended use of the CENTROPIX product, self-made modifications and repairs, misuse, normal wear and tear, catastrophes or force majeure, any claim under warranty is void. The warranty is only guaranteed when using the **KLOUD+** up to a maximum body weight of 120 kg (265 pounds) and when used on a flat, hard surface. Please place the **KLOUD+** on your body if these conditions cannot be met.

Warranty Service

If CENTROPIX identifies a manufacturing defect, the user is entitled to a repair or replacement of a unit under this warranty, provided the damaged product was sent to CENTROPIX within the warranty period. The warranty period will not be extended by any warranty service claimed.

Liability

CENTROPIX is not liable for direct and indirect side effects or consequential and incidental damages, even if these are due to legal, illegal or other actions.

Warranty processing

To exercise your warranty rights in the event of a defect, please contact our customer service department by e-mail at rma@centropix.com or by phone at +423 220 29 59, and let us know your concerns.

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